




MIPROLAB

PERMEABILITY DEVICE

The device is used in soil water permeability tests.
It works with 10 samples at the same time.
It is a closed system in terms of water saving.
It circulates continuously with the help of a pump.
The pump works for 48 hours in water without switching off.
Water tank and water filter.
Design suitable for 53 mm and 60 mm sampling rings.
10 sample rings.
Permeability container made of hard plastic and 53-60 mm diameter seamless stainless steel sampling cylinder made of seamless stainless steel.
Electrostatic powder-coated housing.
220 V, 50 Hz.

**Production of Custom Sizes and Features*

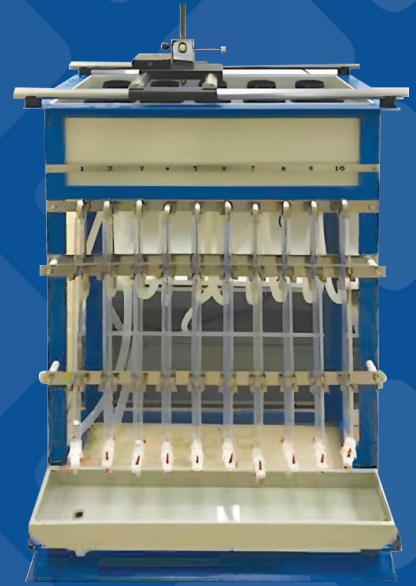

PRESSURE FILTERING DEVICE

It provides fast, reliable and easy extraction of saturation mud extract by air pressure.
It is easy and practical to use and each unit in the system is resistant to 350 kg pressure.
The sealing of the system is provided by special design.



MIPROLAB
MIPROLAB LABORATORY EQUIPMENT INC.

AĞAÇ İŞLERİ SANAYİ SİTESİ 1362. CAD. NO:41 ANKARA/TÜRKİYE

Phone: +90 312 395 75 55 | Whatsapp: +90 542 395 75 56

E-mail: info@miprolab.com.tr

